



ШЛИФОВАЛЬНО-ПОЛИРОВАЛЬНАЯ МАШИНА

QPOL PRO

FLAGSHIP GRINDING AND POLISHING MACHINE - POWER MEETS PRECISION

The Qpol PRO is a fully automatic grinding and polishing machine that operates with both single and central pressure. It was developed for continuous industrial use with maximum operating capacity and impresses with its unique features.

ПРЕИМУЩЕСТВА ОБОРУДОВАНИЯ

- | **Maximum reproducibility:** Material removal measurement with a precision of up to 10 µm, electronic force control in closed-loop, and a sturdy aluminum housing
- | **High performance:** Powerful, maintenance-free direct drives (patented) and sample holders up to 204 mm
- | **Intuitive operation:** Patented sample holder mounting, well-designed software concept, and one-button operation



[Смотреть видео](#)

ШЛИФОВАЛЬНО-ПОЛИРОВАЛЬНАЯ МАШИНА QPOL PRO


ОСОБЕННОСТИ



DURABLE AND RELIABLE

Machine body made of a sturdy welded aluminum construction with a high-quality, durable powder coating

- | Extremely smooth operation thanks to a low-vibration design
- | Maintenance-free direct drives for continuous operation with optimal performance, even for larger samples
- | The impact-resistant plastic tub offers impressive durability, corrosion resistance, and easy cleaning

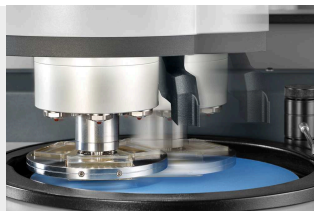

GRINDING/POLISHING HEAD WITH STEP MOTOR

Thanks to its uniform oscillating motion, the grinding/polishing head can utilize the entire working disc

- | Longer service life of consumables
- | Intuitive adjustment of the oscillating motion and ability to memorize previously set endpoints
- | Semi-automatic sample holder mounting (can also be operated with one hand)



PRECISE MATERIAL REMOVAL MEASUREMENT

Thanks to its robust design and force sensors, automatic target preparations can be performed


- | Accuracy of up to 10 μm
- | Closed-loop electronic force control
- | Ideal for electronics and thin sections



ERGONOMICALLY SHAPED SPLASH GUARD

User-friendly design and the ability to quickly replace consumables

- | The ergonomically shaped splash guard significantly reduces splashing and can be used as a comfortable hand rest during longer procedures
- | Grip recesses allow for the replacement of grinding/ polishing media even when the ring is in place

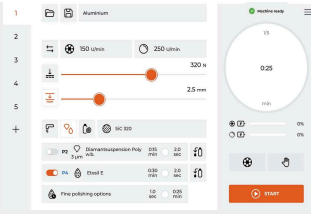

SIMPLE CLEANING SYSTEM

Single-Disc-Design:
Easy removal of all components and easy access

- | The splash guard ring and working wheel can be removed with minimal effort, thereby reducing setup and cleaning times
- | Working disc cooling for long polishing cycles and reduced cleaning effort thanks to basin rinsing
- | The tub insert (plastic inlay) is easy to replace or clean to remove abrasion and residue from the Suspension

INTUITIVE AND EASY TO USE

Ergonomically integrated 10-inch industrial touchscreen with QATM user software

- | Clear menu structure and easy configuration and saving of all process parameters
- | Progress indicator for the current procedure via a status bar and an audible signal at the end of the grinding/polishing process
- | Level indicator for suspensions
- | Pre-installed preparation methods for various materials and applications

OPTION: QDOSER DOSING SYSTEMS

Can be optionally equipped with software-controlled dosing systems for high process flexibility


- | Modular Dosing System: Qdoser GO and PLUS can be combined for up to 8 suspensions
- | Qdoser GO: Fully automatic, software-controlled dispensing system with 4 easily removable and refillable dispensing bottles for reproducible results (magnetic stirring function available as an option)
- | Qdoser PLUS: Also allows for fine polishing (optional upgrade to a 1-liter bottle for fine polishing suspension)

**PRECISE MATERIAL REMOVAL
MEASUREMENT, HIGH-PERFORMANCE
MOTORS AND LARGE SAMPLES**

QPOL PRO - EFFICIENT PREPARATIONS OF THE HIGHEST QUALITY

The Qpol PRO delivers 5–80 N of contact force in single-sample mode and 20–500 N in central pressure mode. Thanks to its powerful direct drive, up to 5 specimens with a diameter of 70 mm can be prepared at variable speeds in synchronous or counter-rotating modes

- | Easy sample holder change, even with just one hand (patented)
- | Reproducible results thanks to stepper motors and material removal measurement
- | Fully controllable via the software; no more swiveling required

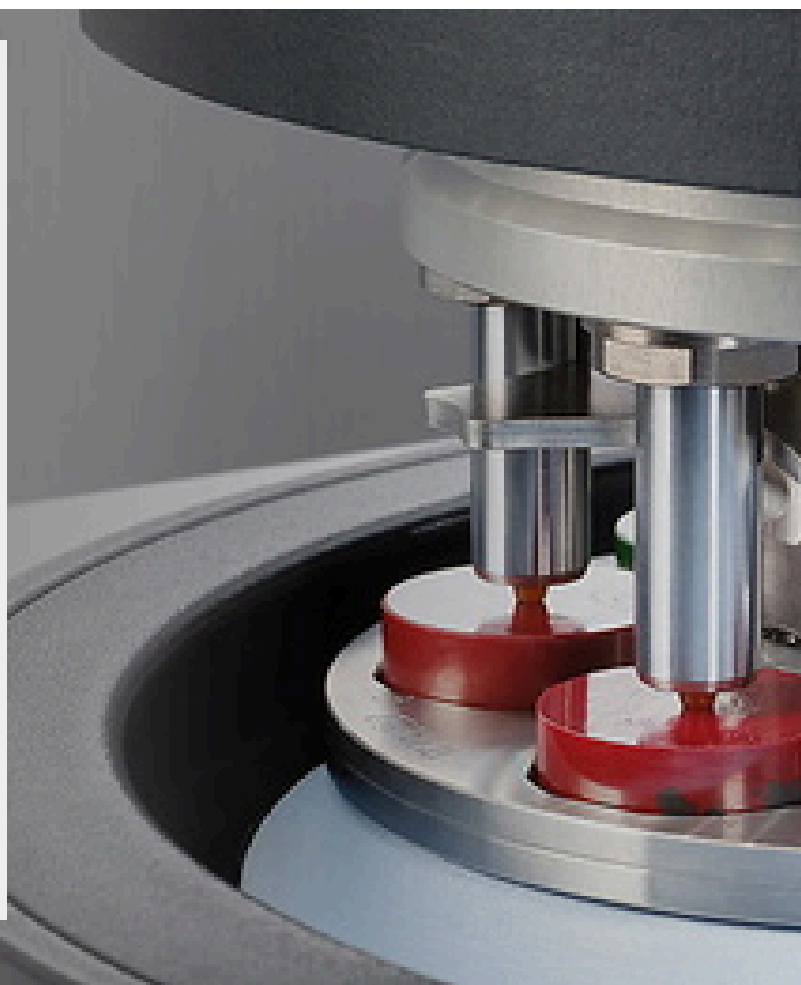

ШЛИФОВАЛЬНО-ПОЛИРОВАЛЬНАЯ МАШИНА QPOL PRO SUPERIORITY IN DETAIL

1. Highest precision

- | Direct force measurement (patented)
- | Material removal measurement system

2. Process-optimized polishing head

- | Intelligent sample holder clamping system (patented)
- | Efficient pendulum motion

3. Robust design

- | Low-vibration cast aluminum components
- | Maintenance-free & high-performance direct drives

4. Intuitive process control

- | 10" industrial touch panel
- | Pre-installed method database
- | One-button operation

5. Modular dosing system

- | Magnetic stirring function
- | Fill level indicator



HIGH PROCESS FLEXIBILITY QDOSER DOSING SYSTEMS

The fully automatic, software-controlled Qdoser dosing system speeds up the grinding and polishing processes in the materialographic laboratory and is compatible with the QPol M, GO, and PLUS series, as well as the QPol PRO. Four screw-top 500-ml glass bottles can be filled with QPREP diamond suspensions, QPREP lubricants (diamond lubricants), and/or QPREP fine-polishing suspensions through a separate opening in the lid. The suspension storage bottles are easily removable, making them simple to clean. A backflush function keeps the connecting hoses clean.






OPTIMIZED FILLING AND DOSING

The suspension storage bottles can be removed at any time and refilled in a way that is optimized for application, e.g., with:

- | Qdoser PLUS: 3x diamond suspensions/ lubricant and 1x fine polishing suspension (optional upgrade to a 1-liter bottle for fine polishing suspension)
- | Qdoser GO: 4x Dia-COMplete diamond suspensions

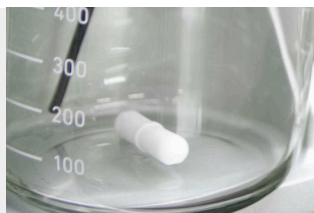

OPTION: MAGNETIC STIRRER

A magnetic stirrer is available as an option for stirring the polishing agents to achieve excellent preparation results

- | The integrated magnetic stirrer prevents the suspension from settling
- | Optimal distribution of polishing particles through automatic stirring of the suspensions



DOSAGE

High flexibility and ease of use thanks to simple polishing agent dispensing

- | Easy installation of the dispensing nozzles into the designat-ed mount in the grinding and polishing head
- | Fully automatic, software-controlled dosing via the industrial touchscreen on Qpol machines
- | Dosing interval can be adjusted as needed
- | Including backflush function for easy cleaning



QDOSER ONE

Dropper bottle for easy and continuous dosing

- | Easy to use – Standard QATM suspension and lubricant bottles can be quickly and easily screwed onto the dispenser
- | Continuously variable dosing rate – The dosing rate is precisely controlled using an adjustment screw
- | Can be used with a stand or mounted on polishing head


PERFECT GRINDING AND POLISHING – CONSUMABLES BY QATM

Grinding and polishing is more reliable, safe and efficient when processed with the right ingredients: knowledge, care and the right choice of consumables and instruments.



ШЛИФОВАЛЬНО-ПОЛИРОВАЛЬНАЯ МАШИНА QPOL PRO

ХАРАКТЕРИСТИКИ

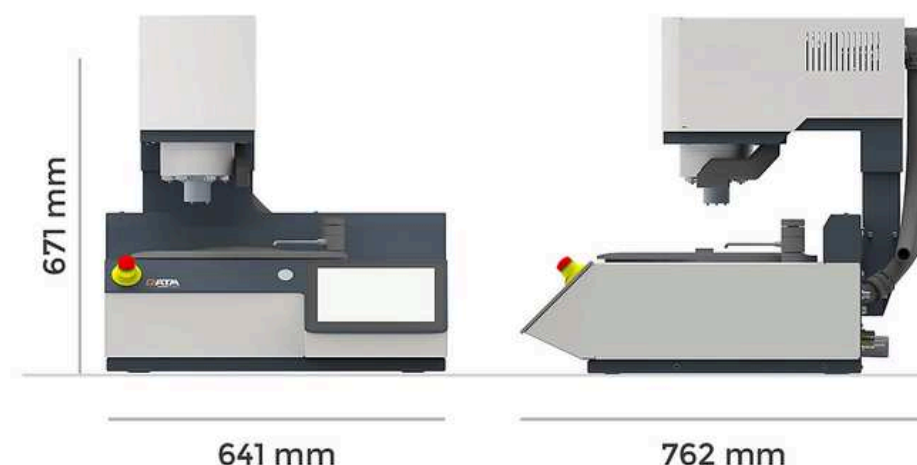

Рабочие диски

ø 250/300 mm
(recommended size: ø 300 mm)

Держатель образцов

Single pressure: ø 160 mm
Central pressure: ø 204 mm

Number of samples

Single pressure: 1 - 6 samples, ø 50 mm
Central pressure: 1 - 5 samples, ø 70 mm

Contact force

Single pressure: variable, 5 - 80 N
Central pressure: variable, 20 - 350 N

Speed (grinding machine)

50 - 600 rpm, continuously variable

Rotational speed (grinding head)

30 - 200 min⁻¹

Подключение воды	1x Fresh water G ½" max. 6 bar
Подключение сжатого воздуха	6-8 bar, coupling connector NG 8
Dosing device (optional)	Qdoser PLUS, Qdoser GO
Температура окружающей среды	10 - 40 °C
Мощность привода (шлифовальный круг)	Working wheel: 1.1 kW (S1) Grinding head: 0.37 kW (S1)
Connected load	1,4 кВА
Power connection	220 - 240 V, 50/60 Hz (1 Ph/N/PE) 100 - 120 V, 50/60 Hz (1 Ph/N/PE)
Внешние интерфейсы	USB, сеть
В x Н x Т	641 x 671 x 762 mm
Вес (в зависимости от комплектации)	approx. 95 kg

www.qatm.com/qpol-pro

ИНФОРМАЦИЯ ДЛЯ ЗАКАЗА

БАЗОВОЕ ОБОРУДОВАНИЕ

МАШИНА (НЕОБХОДИМА - 1 ШТ.)

M5633000		Qpol PRO 200-240 V 50/60 Hz
----------	---	-----------------------------

M5633005		Qpol PRO 100-120 V 50/60 Hz
----------	--	-----------------------------

КОМПЛЕКТ ДЛЯ ПОДКЛЮЧЕНИЯ

Z5631604		Connection set for Qpol PRO
----------	--	-----------------------------

РАБОЧИЙ ДИСК (1 ШТ. НЕОБХОДИМА)

Z5633002		Диск рабочий Ø 250 мм
----------	--	-----------------------

Z5633003		Диск рабочий Ø 300 мм
----------	--	-----------------------

PROTECTIVE COVER (TO BE ORDERED WITH THE MACHINE)

A5633001		Preparation protective cover with safety switch
----------	--	---

Z5633001		Protective cover
----------	--	------------------

ОПЦИИ (МОГУТ ЗАКАЗЫВАТЬСЯ ОТДЕЛЬНО)

		Magnetic foil anisotropic
95016346		Magnetic foil anisotropic Ø 250 x 1.0 mm PSA

95012161		Magnetic foil anisotropic Ø 300 x 1.0 mm PSA
----------	---	--

SPLASHING GUARD (1 ITEM REQUIRED)

Z5631602		Splashing guard Ø 200/250 mm
----------	--	------------------------------

Z5631607		Splashing guard Ø 300 mm
----------	--	--------------------------

CLAMPING AND PROTECTION RING

Z5631004 Кольцо прижимное Ø 250 мм

Z5631009 Кольцо прижимное Ø 300 мм

Z5631022 Защитное кольцо Ø 250 мм

Z5631023 Protection ring Ø 300 mm

GRINDING BOWL INLAY

95018173 Qpol grinding bowl inlay (5-pack)

ДЕРЖАТЕЛЬ ОБРАЗЦОВ (НЕОБХОДИМ - 1ШТ.)

Z5446020 Sample holder L-160-6-Ø-25 single pressure

Z5446022 Sample holder L-160-6-Ø-30 single pressure

Z5446025 Sample holder L-160-6-Ø-40 single pressure

Z5446026 Sample holder L-160-6-Ø-50 single pressure

Z5446021 Sample holder L-160-6-Ø-1" single pressure

Z5446023 Sample holder L-160-6-Ø-1 1/4" single pressure

Z5446024 Sample holder L-160-6-Ø-1 1/2" single pressure

Z5446027 Sample holder L-160-6-Ø-2" single pressure

Z5446028 Sample holder L-160-3-RE-55x30 single pressure

Z5446029 Sample holder L-160-3-RE-60x40 single pressure

Z5446030 Sample holder L-160-12-Ø-25 central pressure

Z5446031 Sample holder L-160-8-Ø-30 central pressure

Z5446034 Sample holder L-160-6-Ø-40 central pressure

Z5446035 Sample holder L-160-5-Ø-50 central pressure

Z5446042 Sample holder L-185-3-Ø-60 central pressure

Z5446044 Sample holder L-204-5-Ø-70 central pressure

Z5446030 Sample holder L-160-12-Ø-1" central pressure

Z5446032 Sample holder L-160-8-Ø-1 1/4" central pressure

Z5446033 Sample holder L-160-6-Ø-1 1/2" central pressure

Z5446035 Sample holder L-160-5-Ø-2" central pressure

Z5446036 Sample holder L-160-3-RE-55x30 central pressure

Z5446037 Sample holder L-160-3-RE-60x40 central pressure

Z5446043	Sample holder L-185-3-RE-70x40 central pressure
Z5446045	Sample holder L-204-3-RE-100x50 central pressure
Z5446038	Sample holder L-160-6-Ø-20-40 central pressure

DISTANCE RING

05446012	Distance ring-50-25
05446015	Distance ring-2"-1"
05446011	Distance ring-50-30
05446014	Distance ring-2"-1 1/4"
05446010	Distance ring-50-40
05446013	Distance ring-2"-1 1/2"

LEVELLING PLATE

05400114	 Levelling plate Ø 160 mm 2 mm
05400246	Levelling plate Ø 160 mm 4 mm

ДОЗАТОР

M6300110	Qdoser PLUS
M6300115	Qdoser PLUS with magnetic stirrer
Z6300070	1 l extension final polish suspension
M6300100	Qdoser GO
M6300105	Qdoser GO with magnetic stirrer
Z5631015	 Qdoser ONE
Z5631014	 Holder for Qdoser ONE
Z5631016	 Подставка для Qdoser ONE

SETTLING TANK / RECIRCULATION UNIT

A5810261		Settling tank 45 L insertable 800
A5810237		Settling tank 45 L insertable 900
A5800496		Settling tank 45 L movable
A5633000		Preparation for recirculation unit for Qpol PRO
Z5800012		Recirculation tank 45 L movable
Z5800013		Recirculation tank 45 L insertable
Z5800014		Pump module recirculation cooling
Z5631605		Qpol control unit for recirculation unit
Z5631606		Qpol shift valve for recirculation unit